
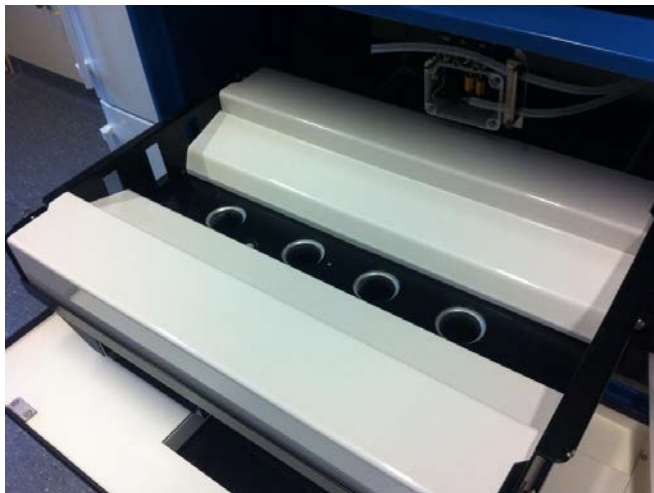


CELL-tainer® expansion-set CBVRC2

The expansion-set makes it possible to create a minimal working volume of 150mL in an easy and reliable way. This is done by placing plastic elements underneath the bag before the bag is placed into the CELL-tainer® tray. As the bag is fully flexible, the shape of the bag when inflating follows the shape of the plastic elements. The four different trays ensure a high flexibility in volume.

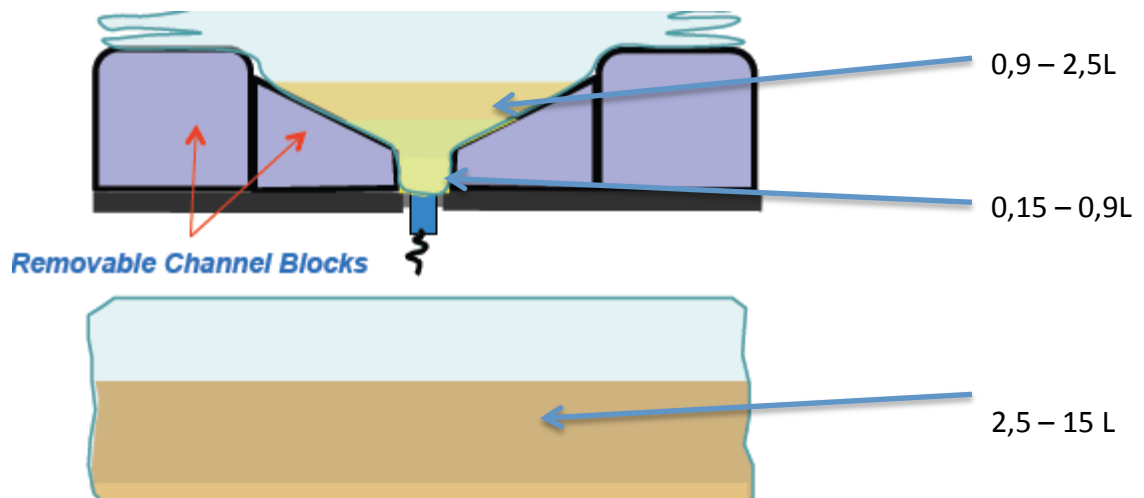

Expansion from 150ml – 15L (or with a 3D square bag up to 25L) can be done in **one and the same bag** saving equipment, bags and handling and reducing risk of contamination. The Expansion set consists of 4 removable plastic elements to be placed underneath the single-use bags.